
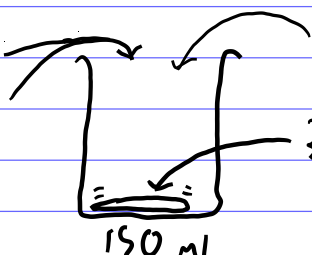
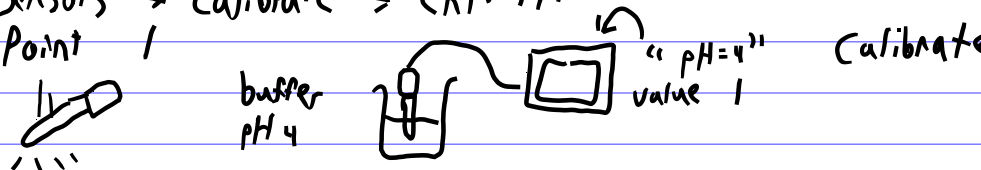
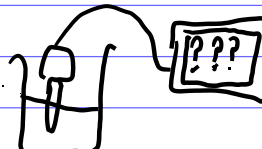


Chem 120 Acid-Base Lab Flowchart


- Part 1
1. 0.35g NaH_2PO_4 & 0.4g Na_2HPO_4
 2. 50 mL DI H_2O
 3. Stir w/ stir bar
- Buffer
- 

- Part 2
- ## Calibration
1. Connect pH meter to Vernier
 2. Add sensor
 3. Sensors → Calibrate → CH1: pH
 4. Point 1
- 
5. Repeat w/ a pH 7 buffer
6. Measure pH of Phosphate buffer
- 

- Part 3
1. add 20 drops of 0.04% bromothymol blue sodium salt indicator solution to buffer & stir

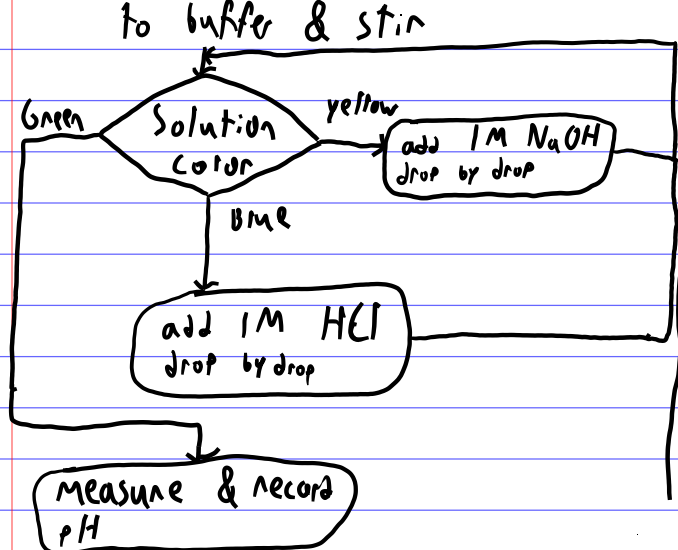

Part 4

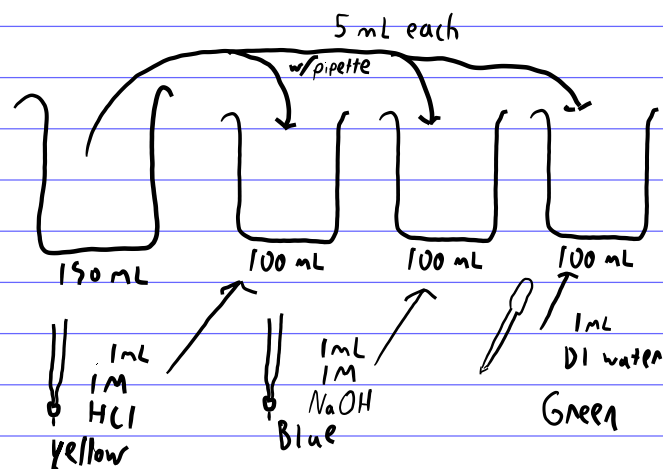

- Part 5
1. Connect Spectrometer
 2. Start Calibration
 3. Insert DI cuvette
 4. Insert yellow cuvette, & measure 453 nm & 616 nm
 5. Repeat for blue & green
 6. Send data

